


# MEDICRAFT BPAP



**MODES:** CPAP, AUTOCPAP, BiLEVEL S  
**PRESSURE RANGE:** 4-25 MBAR  
**MASKS:** Nazal and Oranazal (Optional)  
**ACCESSORIES:** Li-ion Battery, Humidifier, Lighter Adapter  
**OPERATING PRINCIPLE:** Pressure Targeted  
**ADDITIONAL INDICATORS:**  
Tidal Volume Indicator, Frequency, I / E  
**PROGRAMS:** MAP v2  
**Adapter:** 60W

**COLOR:** White  
**WEIGHT:** 1.3 kg  
**DIMENSIONS:** 17 x 23x 23,6 cm  
**SOUND LEVEL:** 28 db  
**ENERGY CONSUMPTION:** 24W  
**POWER SOURCE:** 100-240VAC 50-60HZ 45W  
**STORAGE TEMPERATURE:** -25°C – +75°C  
**OPERATING TEMPERATURE:** 5°C – +40°C  
**HUMIDITY:** 5-95%  
**WARRANTY:** 2 YEARS

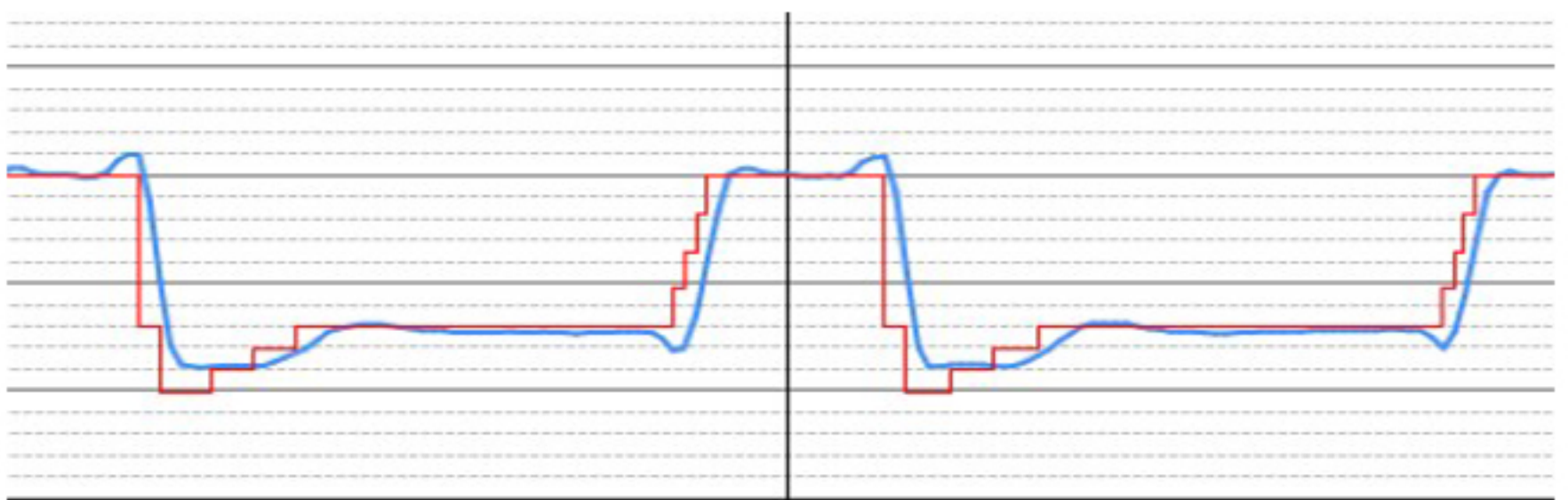

• Reference Pressure

11:35:30

## SLOPE AND PRECISION

The Medcraft slope setting determines the transition time between inspiration and expiration. 4-stage inspiration and 8-stage expiration slope adjustment allows a suitable time to be adjusted. Medcraft devices detect the triggers in the patient's breathing and exhalation passes and apply the inspiratory and expiratory pressures at the right time. These triggers vary according to the lung capacity of the patients.



## BILEVEL FLEXIBILITY FEATURE

There is three stage for flexibility feature. This feature consists of the most important factors that increase the comfort of adaption and therapy.



## 4-25 MBAR OPERATING RANGE

The BiLEVEL device applies inspiratory and expiratory pressures in the range of 4-25 mbar separately. Thus, while the patient takes the pressure patient needs during inspiration, patient can easily breathe with the pressure applied during the expiration without closing the airway.